



# Zen-O™

## Portable Oxygen Concentrator

The Zen-O™ from GCE is an all-in-one portable oxygen solution that caters for a wide range of oxygen patients whilst maintaining a high degree of portability.

It has two oxygen flow modes, providing up to a maximum of 2 litres per minute of continuous flow and 6 pulse flow settings. Weighing from just 10.8lbs (4.66kg), the Zen-O™ is one of very few devices offering a continuous flow option where it's possible that the unit can be worn over the shoulder. The concentrator also comes with a wheeled cart and accessory bag as standard.



Pulse Flow

**1 to 6**  
Pulse flow  
settings



Continuous

**0.5 to 2**  
Litres Per  
Minute



Weight

From  
**10.8lbs**  
(4.6kgs)



Battery life

Up to  
**8**  
Hours  
on setting 2  
pulse flow\*

*intermedical*  
HOMECARE

[www.portableoxygen.co.uk](http://www.portableoxygen.co.uk) | 0800 028 2194

Trusted in Respiratory Health

# Lightweight. Stylish Design. Up to 8 hours battery life.

Zen-O™  
GCE


Its stylish design is coupled with a simple to use interface and intuitive button operation all on a backlit LCD display. Using advanced patented technology, Zen-O™ automatically responds to an increase of the user's breathing rate, increasing the amount of Oxygen delivered to ensure the user remains saturated.

Zen-O™ can be purchased with one or two 12-cell Lithium batteries providing up to 8 hours of use on single charge depending on usage conditions, flow mode and setting. The battery compartment can house two batteries at the same time.



## Single Battery Package

- GCE Zen-O main unit
- 12 cell battery pack
- Pull cart with handle
- Carry case
- AC & DC power supply units
- User manual
- Quick start guide
- Nasal cannula

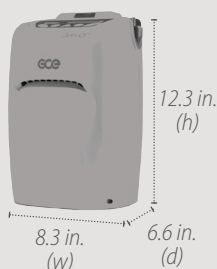

## Double Battery Package


- As above including:
- Additional battery pack


31.3 (h) x 21.2 (w) x 16.8 (d)  
in cm

*\*Run-time on pulse flow mode  
dependent on breathing rate*

Approximate battery run-time  
(in hours:minutes)

	1 x battery pack	2 x battery packs
<b>Continuous Flow Mode</b>		
0.5lpm	3:00	6:00
1lpm	1:45	3:30
1.5lpm	1:15	2:30
2lpm	0:45	1:30
<b>Pulse Flow Mode</b>		
Setting 1	4:00	8:00
Setting 2	4:00	8:00
Setting 3	3:00	6:00
Setting 4	2:15	4:30
Setting 5	2:00	4:00
Setting 6	1:45	3:30

Concentration	90% -3% /+ 6%
Weight	10.8lbs (4.66kg) w/ 1 battery 11.7lbs (5.32kg) w/ 2 batts.
Power	AC: 100-240v DC: 11.5-16v
Noise	38db(A) on setting 2 (pulse)
Batt. charge time	Approx. 4 - 6 hours
Warranty	3 years or 15,000 hours (whichever comes first) 1 year on batteries & accessories
 Airline approved	Yes, subject to airline's terms & conditions

 This device may be suitable for use whilst sleeping. This device may be used with a CPAP or BiPAP device.

**intermedical**  
HOMECARE

0800 028 2194  
[www.portableoxygen.co.uk](http://www.portableoxygen.co.uk)

Intermedical (UK) Limited, Homecare Division  
Unit 6 Mill Hall Business Estate, Aylesford, Kent, ME20 7JZ, United Kingdom



Cert: 87M18032