

## INOGENONE® G4

### Portable Oxygen Concentrator

INOGENONE® G4 breaks the barrier for greater independence and is the perfect solution for those seeking to maintain an active lifestyle whilst on oxygen.

Weighing from as little as 2.8lbs (1.27kg), the INOGENONE® G4 really does “punch above its weight” with an incredibly compact, sleek, and stylish design.

It is the smallest and lightest portable oxygen concentrator in our range offering 3 pulse flow settings.



**Pulse Flow**

**1 to 3**  
Pulse flow  
settings



**Weight**

From  
**2.8lbs**  
(1.27kgs)



**Battery life**

Up to  
**4:30**  
Hours  
on setting 2\*



*intermedical*  
HOMECARE

[www.portableoxygen.co.uk](http://www.portableoxygen.co.uk) | 0800 028 2194

Trusted in Respiratory Health

# The G4 breaks the barrier for greater independence.

The G4 is very quiet and unobtrusive producing noise that is less than a normal indoor conversation. It is simple to operate. The user interface features a white back-lit LCD display that clearly shows the current pulse flow setting as well as the remaining battery life. The pulse flow setting can easily be adjusted using the plus and minus buttons on the control panel.

When using the G4, its intelligent oxygen delivery system works in accordance to the user's respiration rate by automatically adjusting how short or long each burst of oxygen is delivered. During the periods of high exertion, the oxygen is delivered in shorter bursts. Whilst sleeping or during periods when respiration is at a lower rate, the G4 will give longer bursts. This clever system always ensures the user maintains the correct dosage of oxygen as prescribed.

It comes complete with both mains and car power supply units. Also included is a protective carry bag together with a carry strap that can also fix directly onto the machine itself.




◀ The INOGENONE® G4 is available with either a 4 cell battery or an 8 cell battery and slots into the bottom of the machine. A desktop charger is also available to charge batteries independently from the machine.

\* Based on lowest setting

## 4 Cell Battery Package

- INOGENONE® G4 main unit
- 4 cell battery
- AC & DC power supply units
- Carry strap
- Carry bag
- User manual
- Nasal cannula

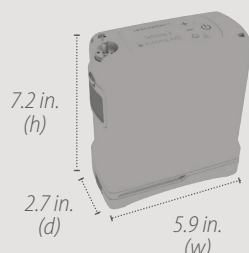

18.2 (h) x 15.0 (w) x 6.8 (d)  
in cm with 4 cell battery installed

## 8 Cell Battery Package

As above but with

- 8 cell battery instead of 4 cell battery

Approximate battery run-time  
(in hours:minutes)

	4 Cell Batt.	8 Cell Batt.
Setting 1	2:40	5:00
Setting 2	2:15	4:30
Setting 3	1:20	2:30

\*Run-time dependent on breathing rate

Concentration	90% -3% /+ 6%
Weight	2.8lbs (1.27kg) with 4 cell battery 3.3lbs (1.49kg) with 8 cell battery
Power	AC: 100-240v DC: 13.5-5.5v

Noise	40db(A) on setting 2
Batt. charge time	3.25hrs with 4 cell battery 6hrs with 8 cell battery
Warranty	3 years. 1 year on batteries, accessories & sieve beds
Airline approved	Yes, subject to airline's terms & conditions
CPAP BiPAP	This device may be suitable for use whilst sleeping. This device cannot be used with a CPAP or BiPAP device.

**intermedical**  
HOMECARE

0800 028 2194  
[www.portableoxygen.co.uk](http://www.portableoxygen.co.uk)

Intermedical (UK) Limited, Homecare Division  
Unit 6 Mill Hall Business Estate, Aylesford, Kent, ME20 7JZ, United Kingdom



Cert: 87M18032