

PHILIPS

RESPIRONICS

SimplyGo

Portable Oxygen Concentrator

Weighing just 10lbs (4.5kgs), the SimplyGo from Philips Respironics easily fulfills the needs of most oxygen users with an all-in-one, continuous and pulse flow oxygen system.

This highly compact concentrator has a pulse flow mode with up to 6 settings and a continuous flow mode that provides a maximum of 2 litres per minute. It also features a unique "sleep mode" function. When activated, pulse flow becomes more sensitive to breath intake and automatically reverts to continuous flow mode if no breath is detected.



Pulse Flow

1 to 6
Pulse flow
settings



Continuous

0.5 to 2
Litres Per
Minute



Weight

10.0lbs
(4.5kgs)



Battery life

**Up to
3**
Hours
on setting 2
pulse flow*



intermedical
HOMECARE

www.portableoxygen.co.uk | 0800 028 2194

Trusted in Respiratory Health

A compact and stylish all-in-one transportable oxygen solution.

SimplyGo is easy to operate with a simple to use control panel and a backlit LCD screen displaying information on current setting, oxygen mode and remaining battery life.

Battery run-time on the SimplyGo is approximately 3 hours on setting 2 pulse flow mode. The optional external battery module allows for another battery to be connected to the machine to extend the overall operation run-time - ideal for long days out or travelling.

SimplyGo comes complete in a stylish brown bag with shoulder strap together with a large wheel cart and accessory bag ideal for keeping all accessories and cables together.




Easy to use control panel with backlit LCD display.



▲ The SimplyGo comes complete with a trendy carry and accessory bag to keep all power cables and other accessories neatly together. A large wheel cart is also included for easier transportation.

Standard Package

- SimplyGo main unit
- Battery pack
- Cart with telescopic handle
- AC mains power supply unit
- DC power supply unit
- Carry case (brown)
- Accessory bag (brown)
- User manual
- Nasal cannula

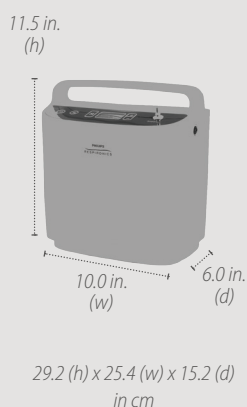

Approximate battery run-time (in hours:minutes)

Continuous Flow Mode		Pulse Flow Mode	
0.5 lpm	2:50	Setting 1	3:20
1 lpm	1:35	Setting 2	3:00
2 lpm	0:50	Setting 3	2:10
		Setting 4	1:40
		Setting 5	1:15
		Setting 6	1:15

*Run-time dependent on breathing rate on pulse flow mode

Concentration	90% -3% /+ 6%
Weight	10.0lbs (4.5kgs) with battery installed
Power	AC: 100-240v DC: 19v
Noise	43db(A) on setting 2 (pulse)
Batt. charge time	Approx. 2 - 5 hours
Warranty	2 years 90 days on sieve beds, battery & accessories
✈️ Airline approved	Yes, subject to airline's terms & conditions



This device may be suitable for use whilst sleeping. This device may be used with a CPAP or BiPAP device.

intermedical
HOMECARE

0800 028 2194
www.portableoxygen.co.uk

Intermedical (UK) Limited, Homecare Division
Unit 6 Mill Hall Business Estate, Aylesford, Kent, ME20 7JZ, United Kingdom



Cert: 87M18032

Trusted in **Respiratory Health**