

NEW!

INOGENONE® G5

Portable Oxygen Concentrator

The INOGENONE® G5 is a remarkable feat of medical engineering, offering users an exciting technological innovation in portable oxygen therapy with a high capacity output, in a truly compact and stylish design.

In a first of its kind, the INOGENONE® G5 meets the needs of a wide spectrum of oxygen users with 6 pulse flow settings in an incredibly small and lightweight solution weighing from just 4.7lbs (2.1kg).

The INOGENONE® G5 can be fitted with an 8-cell or a 16-cell battery providing up to 6 ½ or 13 hours of operation time respectively, making it perfect for travel as well as everyday use.



Pulse Flow

1 to 6
Pulse flow
settings



Weight

From
4.7lbs
(2.1kgs)



Battery life

Up to
10
Hours
on setting 2*



intermedical
H O M E C A R E

www.portableoxygen.co.uk | 0800 028 2194

Trusted in Respiratory Health

Powerfully Portable.

In operation, the INOGENONE® G5 is very quiet and no more obtrusive than the sound of a normal indoor conversation. The provided black protective bag is of a stylish and durable material made especially for keeping the INOGENONE® G5 as discreet as possible.

The LCD display and backlit buttons on the user interface are simplistic and self-explanatory making it very easy to use.

Designed for 24 hour use, the INOGENONE® G5 is built with INOGEN's intelligent pulse oxygen delivery system. This technology works in conjunction with the user's respiration rate that automatically adjusts the amount of oxygen that is delivered per breath. In periods of rest or sleep, the INOGENONE® G5 will cleverly deliver longer bursts of oxygen in order to retain the user's prescribed flow rate.




◀ The INOGENONE® G5 is available with either an 8-cell battery or a 16-cell battery and slots into the bottom of the machine. A desktop charger is also available to charge batteries independently from the machine.

* Based on lowest setting

8 Cell Battery Package

- INOGENONE® G5 main unit
- 8 cell battery
- AC & DC power supply units
- Carry strap
- Carry bag
- User manual
- Nasal cannula

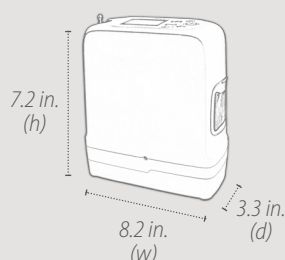

18.2 (h) x 20.7 (w) x 8.4 (d)
in cm with 8 cell battery installed

16 Cell Battery Package

- As above but with
- 16 cell battery instead of 8 cell battery

Approximate battery run-time
(in hours:minutes)

	8 Cell Batt.	16 Cell Batt.
Setting 1	6:30	13:00
Setting 2	5:00	10:00
Setting 3	3:30	7:00
Setting 4	2:30	5:00
Setting 5	1:40	3:20
Setting 6	1:15	2:30

*Run-time dependent on breathing rate

Concentration	90% -3% /+ 6%
Weight	4.7lbs (2.1kg) with 8 cell battery 5.3lbs (2.4kg) with 16 cell battery

Power	AC: 100-240v DC: 13.5-15v
Noise	37db(A) on setting 2
Batt. charge time	3hrs with 8 cell battery 6hrs with 16 cell battery
Warranty	3 years. 1 year on batteries, accessories & sieve beds
Airline approved	Yes, subject to airline's terms & conditions



This device may be suitable for use whilst sleeping. This device cannot be used with a CPAP or BiPAP device.

intermedical
HOMECARE

0800 028 2194
www.portableoxygen.co.uk

Intermedical (UK) Limited, Homecare Division
Unit 6 Mill Hall Business Estate, Aylesford, Kent, ME20 7JZ, United Kingdom

